

## LITTY 4all

THE LIGHT & STURDY  
CHILDREN'S WHEELCHAIR THAT  
GROWS WITH ITS ADVENTURER.

The LITTY 4all **children's wheelchair** is specially designed for little adventurers. With its optimal frame geometry and versatile adjustment options, it adapts perfectly to the needs of children. Its low weight from 6.7 kg and the particularly low rolling resistance make it the ideal companion even for children with little arm strength. With the LITTY 4all, 9 frame parameters can be selected to customise the wheelchair, and the seat width is possible from as little as 18 cm.




### #GrowAlongFeature

The LITTY 4all grows with the child thanks to interchangeable assemblies and thus remains a reliable companion for years.



### #PerfectSittingPosition

The frame structure, without restricting front frame tubes, allows for seat positioning with maximum abduction.

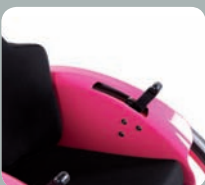

### #Safety

The operation of the knee lever brake is ergonomically mounted in the clothing protection and the brake bolt is outside the grip area to prevent fingers from getting trapped.



### #Support

Various push handles and bars offer the possibility of assistance when moving forward, while the companion brakes allow speed control.





## PERFECTLY COORDINATED AND ALWAYS READY FOR THE NEXT ADVENTURE.

The **LITTY 4all** is individually manufactured according to the needs of the child – from the frame dimensions to the features. There are a variety of options to choose from, including footrests, brakes, push handles, back and seat systems, and drive and caster wheels.

### UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



This “foot support that folds up to the rear” makes it easier to transfer in and out of the wheelchair for children who are still able to stand.



The therapy table offers a practical work surface, can be removed without tools and can be infinitely adjusted in height and angle for optimal adjustment.



With the **outdoor wheel attachment**, the LITTY 4all is suitable for off-road use on pavements, gravel and forest paths. The caster wheels are lifted off the ground and the 12" wheel overcomes bumps effortlessly. Thanks to the quick-release clamp, the wheel attachment can be attached easily and without tools.



## TECHNICAL SPECIFICATIONS

Seat width:	18 - 34 cm (in 2 cm steps)
Seat depth:	20 - 35 cm (in 5 cm steps)
Back height:	20 - 35 cm (in 5 cm steps)
Lower leg length:	20 - 44 cm (in 1 cm steps)
Rear seat height:	31 - 41 cm (in 1 cm steps)
Front seat height:	34 - 44 cm (in 1 cm steps)
Tipping point:	5 - 9 cm (in 1 cm steps)
Total product weight:	Starting at 6.7 kg
Max. standard load weight:	Up to 50 kg

Data subject to design changes.

## FOLLOW US ON SOCIAL MEDIA



**PROACTIV®**

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11  
72359 Dotternhausen  
info@proactiv-gmbh.de

[www.proactiv-gmbh.com](http://www.proactiv-gmbh.com)

Rev. 09/2024\_01

Subject to technical modifications, printing errors & colour deviations.